

Saving a life deserves the best protection



Defib Store 4000



**FIRE
RETARDANT**



**FULLY
RECYCLABLE**



**WATER
RESISTANT
IP66**



**ROBUST
DESIGN**

Revolutionary polycarbonate defibrillator cabinet from Ce-Tek, the industry leading AED cabinet manufacturer. Uncompromising on quality. Uncompromising on safety.

- Highly impact resistant polycarbonate outdoor cabinet
- Fire retardant to UL94 V0 (will not support fire)
- UV stabilised; will not be affected by sunlight
- Double skin clamshell design with full IP66 rating
- Marine grade stainless steel keypad lock
- Thermostatically controlled heater
- Motion sensor LED light
- Lightweight construction, only 10kg
- Fully recyclable material
- 10 Year limited guarantee
- Made in the UK



www.defibstore.co.uk

Defib Store recognise the importance of a high-quality, secure cabinet to store your public access defibrillator

Manufactured from polycarbonate the innovative Defib Store 4000 is the only polycarbonate AED cabinet on the market that can withstand high impact, is fire retardant to UL94 V0 and is UV stabilised against the effect of direct sun. The Defib Store 4000 is the latest innovative cabinet addition to the widely acclaimed Defib Store range of cabinets offering a durable, secure outdoor AED cabinet suitable for all environments.

Defib Store 4000 features include:

Unique double skin construction provides a robust impact-resistant structure plus creates additional insulation to the defibrillator.

Moulded from fire retardant polycarbonate to UL 94 V0, the Defib Store 4000 is self-extinguishing and does not support fire. The polycarbonate is also UV stabilised and will not degrade in direct sunlight.

Fully sealed clamshell design, the sleek cabinet is IP66 EN60529 rated providing a completely dust and water tight AED cabinet.

Thermostatically controlled heater ensures the AED is always within its working temperature range.

Marine grade, stainless steel mechanical keypad lock; tested with salt spray for 1,000 hours, ensures the lock remains fully functional in all weather. Also available unlocked with a plastic winged handle and optional alarm*.


Battery powered, motion sensor LED light to illuminate cabinet on opening. Easily detachable for use as a torch in a rescue.




Defib Store 4000 Ile



Defib Store 4000 Ile

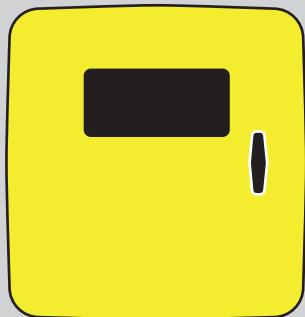

Defib Store 4000 ule

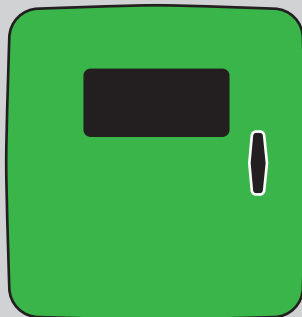

Defib Store 4000 ule

Internal AED hook on back wall accommodates any defibrillator ensuring AED is elevated from base and clearly visible through large viewing window.

Sloping design to allow surface water to drain away from the door.

Available in high visibility yellow or green as standard. Other colours available subject to minimum order quantities.

All Defib Store cabinets include printed graphics on high grade laminated vinyl with the option to include personalised logos/text.

Simple to install with dual exit point for power cable via rear or base of cabinet without onsite drilling or modification.

Cabinets are delivered fully assembled with UK plug, fixing kit, installation and maintenance guide.

10 year limited guarantee.

Defib Store 4000 also available with cooling facility – ask Sales for further information.

Cabinet External Dimensions: 465mm (H) x 435mm (W) x 240mm (D)

Product code	Material	Colour	Door Locking	Heating & Lighting	IP Rating	Viewing Window	Location
DFS4K442LLE-Y	Polycarbonate	Yellow	SS Keypad lock	Yes	IP66	Yes	Indoor/Outdoor
DFS4K442LLE-G	Polycarbonate	Green	SS Keypad lock	Yes	IP66	Yes	Indoor/Outdoor
DFS4K442ULE-Y	Polycarbonate	Yellow	Plastic Wing Handle	Yes	IP66	Yes	Indoor/Outdoor
DFS4K442ULE-G	Polycarbonate	Green	Plastic Wing Handle	Yes	IP66	Yes	Indoor/Outdoor

*Small additional charge



www.defibstore.co.uk

Defib Store Ltd, Units 1-4 Tideswell Business Park, Tideswell, Derbyshire SK17 8NY

Tel: +44 (0) 1298 872186

Fax: +44 (0) 1298 871155

E-mail: sales@defibstore.co.uk